


SOL-CARE™ Safety Needle

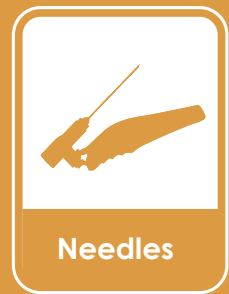

The SOL-CARE™ Safety Needle is designed to prevent needlestick injury.

Its robust safety mechanism is intuitive and compatible with all standard luer lock and slip tip syringes.


Available in a broad range of lengths and gauges.

- Easy, one-handed activation of the protective arm.
- Audible click indicates activation of the safety mechanism.
- Patented locking mechanism ensures that the needle is secure and locked inside the protective arm.
- Hub, protective arm, and packaging are color coded for easy and fast identification of 18G-30G needles (ISO 6009).

Two options for activation of the protective arm:




activation against hard surface

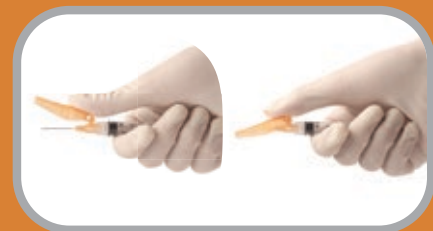


activation with thumb



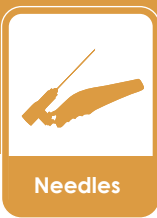

SOL-CARE™ Safety Needle

- In compliance with criteria given in the **EU Directive 2010/32/EU**.
- Protective arm **doesn't obstruct visualization of injection site** (1).
- **Oriented bevel-up needle** allows for low-angle injections (2).
- Hub, protective arm and packaging are **color coded** for easy and fast identification **18G -30G needles** (ISO 6009).
- Available in 8 needle length **from 10 to 50 mm**.
- **Compatible with all standard** luer lock and slip tip **syringes**.

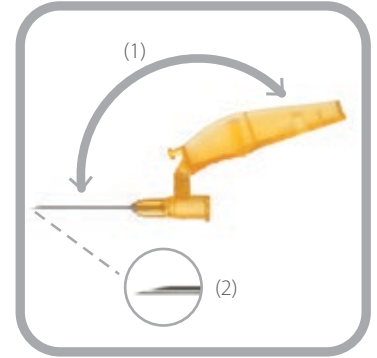

Reference number	Product description	Gauge size	Outer diameter	Length		Indications*				
SN3012	SOL-CARE™ Safety Needle 30G*1/2"	30G	0.3 mm	12.5 mm	1/2"	SC				
SN2712	SOL-CARE™ Safety Needle 27G*1/2"	27G	0.4 mm	12.5 mm	1/2"	SC	ID			
SN2515	SOL-CARE™ Safety Needle 25G*1 1/2"	25G	0.5 mm	40 mm	1 1/2"	SC	ID	IM		IA
SN2510	SOL-CARE™ Safety Needle 25G*1"	25G	0.5 mm	25 mm	1"	SC	ID	IM		IA
SN2558	SOL-CARE™ Safety Needle 25G*5/8"	25G	0.5 mm	16 mm	5/8"	SC	ID	IM		IA
SN2310	SOL-CARE™ Safety Needle 23G*1"	23G	0.6 mm	25 mm	1"			IM	IV	IA
SN2215	SOL-CARE™ Safety Needle 22G*1 1/2"	22G	0.7 mm	40 mm	1 1/2"			IM	IV	IA
SN2210	SOL-CARE™ Safety Needle 22G*1"	22G	0.7 mm	25 mm	1"			IM	IV	IA
SN2115	SOL-CARE™ Safety Needle 21G*1 1/2"	21G	0.8 mm	40 mm	1 1/2"			IM	IV	IA
SN2110	SOL-CARE™ Safety Needle 21G*1"	21G	0.8 mm	25 mm	1"			IM	IV	IA
SN1815	SOL-CARE™ Safety Needle 18G*1 1/2"	18G	1.2 mm	40 mm	1 1/2"				IV	

SC - subcutaneous, ID - intradermal, IM - intramuscular, IV - intravenous, IA - intra-arterial.
 * **Recommended indications for specific measures. Users are explicitly responsible for correct application.**

Sales Unit: 100 pcs/Box; 800 pcs/Case

For more information, please contact:

Sol-Millennium Medical Inc.
 315 Shawnee North Drive, Suite 100,
 Suwanee, GA 30024
 T: 404-973-2200
 e-mail: info@sol-m.com

CORP-SIN-119-0814

sol-m.com

